


KliniMark™ Clip **PRO**

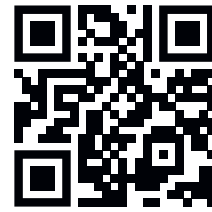

klinimark.com

KliniMark™ Clip PRO Tissue Marker with One-Handed Application System


The KliniMark™ Clip PRO is an innovative marker clip for reliable and permanent soft tissue marking. It has been specifically developed for preoperative use as well as for marking and tracking lesions before or after chemotherapy or radiotherapy.

Based on the proven KliniMark technology, the PRO version offers enhanced ease of use through an ergonomic one-handed application system.

Key Feature

DESIGNED FOR ONE-HANDED USE

Designed for precise placement with one hand, while keeping the other free for ultrasound guidance.
Enables fast, safe and accurate marking — even under challenging conditions.

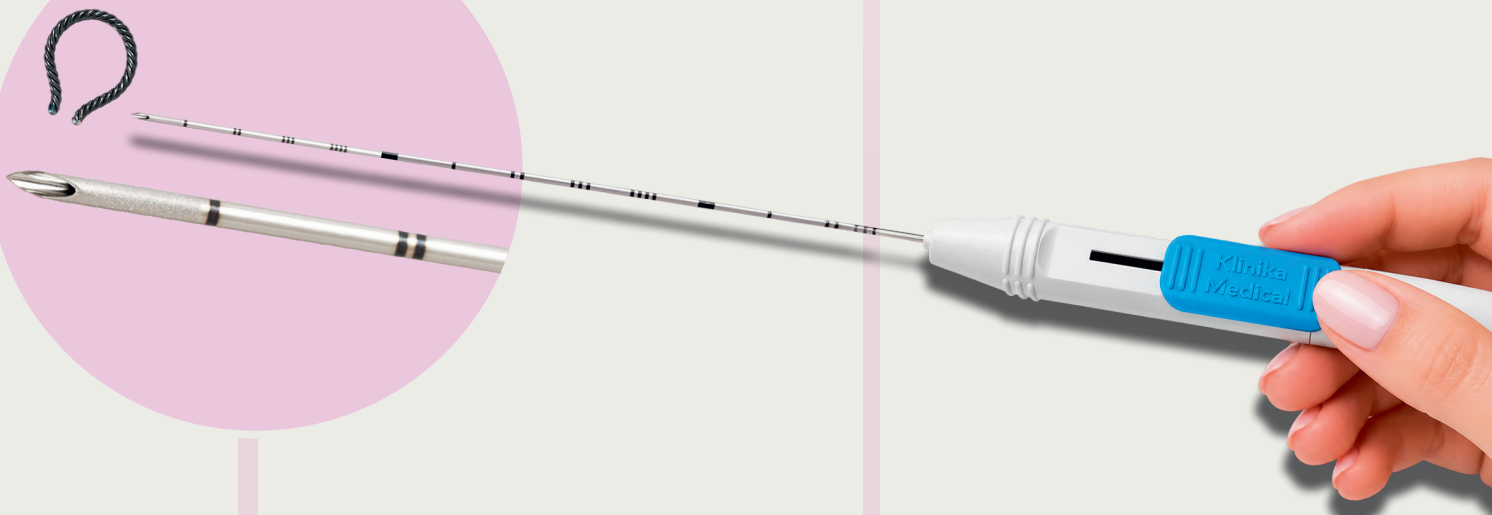

100%
AVAILABILITY



MADE IN
GERMANY



GLOBAL
STANDARD



YOUR LOCAL
SUPPORT



FAMILY-OWNED
COMPANY

➤ Benefits

- **Reliable and permanent tissue marking**
- **Ergonomic one-handed application system for ease of use**

Excellent visibility in all relevant imaging modalities:

- Ultrasound
- Mammography (X-ray)
- MRI
- **Secure fixation in tissue due to the special omega shape**
- **Smooth and controlled placement for high user comfort**
- **Clip is MRI-compatible – clip can remain in the tissue until surgical removal**
- **Biocompatible material (implant-grade nitinol)**

➤ Indications

- Marking of lesions or tumors before or after chemotherapy or radiotherapy
- Marking of biopsy sites
- Identification of resected tumor areas

Preoperative tissue marking to enable precise localization for:


- surgical procedures
- targeted radiotherapy

➤ Application

- KliniMark™ Clip PRO is placed percutaneously via an application needle and deployed in the target tissue.
- Its three-dimensional structure provides excellent visibility in imaging and supports precise therapy planning.

➤ KliniMark™ Clip PRO

Design	KliniMark™ Clip PRO sterile steel cannula with Nitinol marker
Article no.	205880900C
Diameter (gauge)	18
Length of cannula (mm)	155
Cannula material	Steel
Size of clip (mm)	4,9 x 2,1
Diameter of clip (mm)	4,5
Clip material	Nitinol
PU	2x5
Sterilization method	EO/ETO

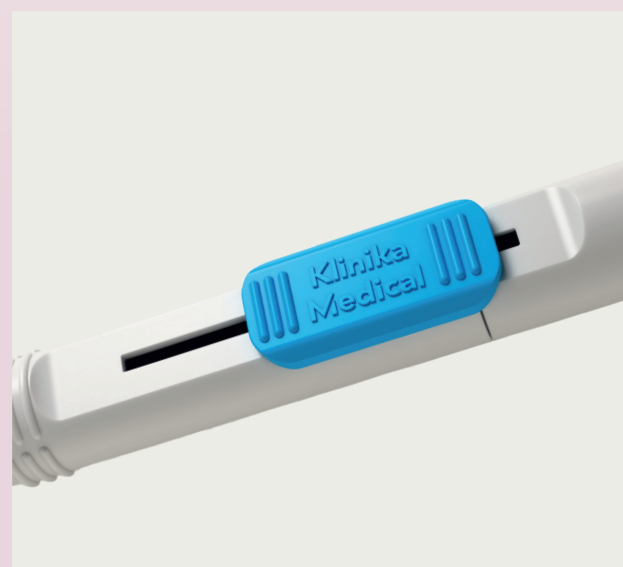

Notes

➤ KliniMark™ Clip PRO

- Sterile single-use product – for one-time use only
- The marker is non-resorbable and remains in the tissue until planned removal
- Contains nickel (nitinol) – use in patients with known nickel allergy only after medical evaluation

